



MASKS TYPE IIR:

<https://drive.google.com/file/d/1PLh8mY982ui3N7P9zjbNh9JfdLHA6ubx/view?usp=sharing>

PRODUCT DOCUMENTS BELOW



DISHANG

Weiwei Dishang Medical Technology Co., Ltd.

Face Mask Non-Sterile Type IIR

Technical Data Sheet

Device variants:

- Type IIR: Disposable Mask with elastic ear loops.
- EN14683:2019



- Shelf life: 2 years
- Storage: Protect from sunlight, dust and humidity
- Usage: Single use
- Sterilization: Non- Sterile
- Applied standards: EN 14683:2019, CE class 1
- Product classification: Medical device acc. to MDD 93/42/EEC.

Measurements:

Item:	Dimensions:	Information:
Mask:	Length	175 +/-
	Width (pleated)	95 +/-
Loop:	Length	175 +/-

Technical Features :

Performance:	Information: Type IIR
Bacterial filtration efficiency (BFE), (%)	≥ 98%
Differential pressure (Pa/cm ²)	< 60
Splash resistance (kPa)	≥ 16,0

Material Data:

Section:	Item:	Information:
Inner layer:	Material	Polypropylene (PP) Non-woven
	Weight/m ²	25 g
	Colour	White
Outer layer:	Material	Polypropylene (PP) Non-woven
	Weight/m ²	25 g
	Colour	Blue
Filter media:	Material	Melt-blown
	Weight/m ²	25 g
	Colour	White
Nose piece:		Pliable encapsulated
Free of:		Glass fibres, latex, additives

IMPORTANT NOTICE:

This guide is only an outline. It should not be used as the only means for selecting protective clothing. Before using any protective clothing, the wearer must read and understand the user instructions for each product.

Specific country legislation must be observed. If in doubt, contact a safety professional. Selection of the most appropriate PPE will depend on the particular situation and should only be made by a competent person knowledgeable of the actual working conditions and the limitations of PPE.

Final determination as to the suitability of these products for a particular situation is the employer's responsibility.

This information is subject to revision at any time. Always read and follow instructions.

LIMITATION OF LIABILITY:

Except as provided above, Dishang Medical shall not be liable or responsible for any loss or damage, whether direct, indirect, incidental, special or consequential arising out of the sale, use or misuse of this product, or the user's inability to use such products.



Dishang Medical, manufacturer of the below referenced product hereby declares that the product:

DS02 Blue Surgical Face Masks Type IIR

Is in conformity with the provisions of European Medical Devices Directive (MDD) 93/42/EEC and, where such is the case, with the national transposing harmonised standard:

EN 14683: 2019 (Type IIR)

As issued by -

SGS-CSTC Standard Technical Services (Shanghai) Co. Ltd
3rd Building No. 889
Yishan Road
Xuhui District
Shanghai
China

Certificate Number: SL52025233855601TX

The Technical Construction file is maintained by Dishang Medical, The Dishang Group
186 West Wenhua Road, Weihai, Shandong, 264209, China.

The product referenced above is identical to the model which is the subject of the EC
Certificate issued by:

Weihai Dishang Medical Technology Co. Ltd
Room 406-409, Block C
No. 213 Torch Road
Torch High-Tech Industrial Development Zone
Weihai
Shandong
China

And further states that it is issued under the sole responsibility of Weihai Dishang Medical
Technology Co Limited.

Completed at: Dishang Medical, The Dishang Group, 186 West Wenhua Road, Weihai,
Shandong, 264209, China.

Eric Wei
Director

Issue Date: 13.05.2020

Package Information of Disposable Medical Surgical Mask

Box size :L19.5cmX 10cm X 8cm

50PCS/inner bag , 1 pre-pack /box , 60 boxes in one ctn



Ctn size :L60cmX 42 cm X 42cm





CERTIFICATION OF REGISTRATION

2020

This certifies that:

**WEIHAI DISHANG MEDICAL TECHNOLOGY CO.,LTD
ROOM406-409 BLOCK C,NO.213 TORCH ROAD,TORCH HIGH-TECH
INDUSTRIAL DEVELOPMENT ZONE, WEIHAI, SHANDONG, 264200, CHINA**

Was registered with US Food and Drug Administration, Center for Devices and Radiological Health, pursuant to the Code of Federal Regulations 21 CFR 807, by GST (Shanghai) Certification Co., Ltd.

Owner/Operator Number:10062870 / Registration Number: 3016680721


Device Listing#:See annex

Expiration Date: December.31,2020

2020 GST will confirm that such registration remains effective upon request and presentation of this certificate until the end of the calendar year stated above, unless said registration is terminated after issuance of this certificate. GST makes no other representations or warranties, nor does this certificate make any representations or warranties to any person or entity other than the named certificate holder, for whose sole benefit it is issued. This certificate does not denote endorsement or approval of the certificate-holder's device or establishment by the U.S. Food and Drug Administration. GST assumes no liability to any person or entity in connection with the foregoing.

Pursuant to 21 CFR 807.39, "Registration of a device establishment or assignment of a registration number does not in any way denote approval of the establishment or its products. Any representation that creates an impression of official approval because of registration or possession of a registration number is misleading and constitutes misbranding." The U.S. Food and Drug Administration does not issue a certificate of registration, nor does the U.S. Food and Drug Administration recognize a certificate of registration, GST is not affiliated with the U.S. Food and Drug Administration.

**for and on behalf of
Greatesting Service Shanghai Co.,Ltd**

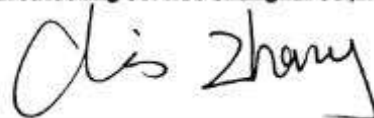

Dated: Aug 3, 2020



CERTIFICATION OF REGISTRATION


2020

Annex to Device Listing#

Proprietary Name	Product Codes	Device Class	Listing Number	Establishment Operations
<i>Disposable Medical Operating Mask</i>	<i>QKR</i>	<i>Not Classified</i>	<i>D384601</i>	<i>Manufacturer</i>

End of the Annex

for and on behalf of
Greatesting Service Shanghai Co.,Ltd

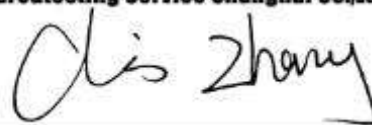

Dated: *Aug 3, 2020*

Establishment Registration & Device Listing

[FDA Home](#) [Medical Devices](#) [Databases](#)

New Search		Back To Search Results
Proprietary Name:	Disposable medical Mask; Disposable Medical Operating Mask; Disposable Protective Mask	
Classification Name:	FACE MASK (EXCEPT N95 RESPIRATOR) FOR GENERAL PUBLIC/HEALTHCARE PERSONNEL PER IIE GUIDANCE	
Product Code:	QKR	
Device Class:	Not Classified	
Registered Establishment Name:	WEIHAI DISHANG MEDICAL TECHNOLOGY CO., LTD	
Registered Establishment Number:	3016990721	
Owner/Operator:	Weihai Dishang Medical Technology Co., Ltd	
Owner/Operator Number:	10062570	
Establishment Operations:	Manufacturer	

Page Last Updated: 08/10/2020

Note: If you need help accessing information in different file formats, see [Instructions for Downloading Viewers and Players](#).Language Assistance Available: [Español](#) | [繁體中文](#) | [Tiếng Việt](#) | [한국어](#) | [Tagalog](#) | [Русский](#) | [العربية](#) | [Kreyòl Ayisyen](#) | [Français](#) | [Polski](#) | [Português](#) | [italiano](#) | [Deutsch](#) | [日本語](#) | [فارسی](#) | [English](#)[Accessibility](#)[Contact FDA](#)[Careers](#)[FDA Basics](#)[FOIA](#)[No FEAR Act](#)[Nondiscrimination](#)[Website Policies](#)

U.S. Food and Drug Administration

10903 New Hampshire Avenue
Silver Spring, MD 20993
Ph. 1-866-INFO-FDA (1-866-463-6332)
[Contact FDA](#)[For Government](#)[For Press](#)[Combination Products](#)
[Advisory Committees](#)
[Science & Research](#)
[Regulatory Information](#)
[Safety](#)
[Emergency Preparedness](#)
[International Programs](#)
[News & Events](#)
[Training and Continuing Education](#)
[Inspections/Compliance](#)
[State & Local Officials](#)
[Consumers](#)
[Industry](#)
[Health Professionals](#)
[FDA Archive](#)

U.S. Department of Health & Human Services

Test Report **SL52025233855601TX** **Date: April 20, 2020** **Page 1 of 3**
 WEIHAI DISHANG MEDICAL TECHNOLOGY CO., LTD
 ROOM 406-409, BLOCK C, NO.213 TORCH ROAD, TORCH HIGH-TECH INDUSTRIAL DEVELOPMENT
 ZONE, WEIHAI, SHANDONG CHINA

The following sample(s) was/were submitted and identified on behalf of the client as:

Sample Description : (A) Disposable medical surgical mask (Claimed Type IIR)
 Sample Color : (A) white
 Style No. : DS01-WHITE/DS02-BLUE
 Proposed Care Instruction : -
 Test Performed : Selected test(s) as requested by applicant
 Sample Receiving Date : Mar 30, 2020 & Apr 02, 2020
 Testing Period : Mar 31, 2020 - Apr 20, 2020
 Test Result(s) : For further details, please refer to the following page(s).

Comment:

Medical Face Masks-Requirements and Test Methods(EN 14683:2019)	(A)
Clause 5.2.2 Bacterial filtration efficiency (BFE)	M
Clause 5.2.3 Breathability	M
Clause 5.2.4 Splash Resistance	M
Clause 5.2.5 Microbial Cleanliness	M

Remark: M=Meet EN 14683:2019 Type IIR requirement

Signed for and on behalf of
 SGS-CSTC Standards Technical Services (Shanghai) Co., Ltd Testing Center

Sara Guo

Sara Guo (Account Executive)



Unless otherwise agreed in writing, this document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Document.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law. Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Attention: To check the authenticity of testing / inspection report & certificate, please contact us at telephone: (86-755) 8307 1443, or email: CN.Doccheck@sgs.com

3rd Building, No. 889, Yishan Road, Xuhui District Shanghai, China 200233
 中国·上海·徐汇区宜山路889号3号楼 邮编: 200233

t (86-21) 61402666 f (86-21) 64958763
 t (86-21) 61402666 f (86-21) 64958763

www.sgsgroup.com.cn
sgs.china@sgs.com

Test Result

Medical Face Masks-Requirements and Test Methods

(EN 14683:2019)

Clause 5.2.2 Bacterial filtration efficiency (BFE)@**

	1#	2#	3#	4#	5#
(BFE), %	99.6	99.9	99.8	99.7	99.6

Remark: Performance Requirement: Type I > 95%, Type II > 98%, Type IIR > 98%

** : The test was carried out by external laboratory assessed as competent

@ : These test methods are not in CNAS accredited scope

Clause 5.2.3 Breathability

(EN 14683 :2019 Annex C, Flow rate 8 l/min)

Sample: A

	1#	2#	3#	4#	5#
Differential pressure ΔP (Pa/cm ²)	32.7	29.4	30.1	28.7	26.5

Remark: Performance Requirement: Type I < 40 Pa/cm², Type II < 40 Pa/cm², Type IIR < 60 Pa/cm²

Clause 5.2.4 Splash Resistance

(ISO 22609 :2004, Pressure 16.0 kPa)

Sample: A

Penetration on inside surface

1#	2#	3#	4#	5#	6#	7#	8#
Pass	Pass	Pass	Pass	Pass	Pass	Pass	Pass
9#	10#	11#	12#	13#	14#	15#	16#
Pass	Pass	Pass	Pass	Pass	Pass	Pass	Pass
17#	18#	19#	20#	21#	22#	23#	24#
Pass	Pass	Pass	Pass	Pass	Pass	Pass	Pass
25#	26#	27#	28#	29#	30#	31#	32#
Pass	Pass	Pass	Pass	Pass	Pass	Pass	Pass

Number of Pass: 32
Overall result: Acceptable

Remark:

- 1) Performance Requirement Type I: N/A, Type II: N/A, Type IIR: ≥ 16.0 kPa
- 2) Distance of the medical face mask target area surface to the tip of cannula is 300 ± 10 mm.
- 3) Condition and Test temperature $(21 \pm 5)^\circ \text{C}$, relative humidity $(85 \pm 10)\%$
- 4) An acceptable quality limit of 4.0% is met for a single sampling plan when 29 or more of the 32 tested specimens show pass results



Unless otherwise agreed in writing, this document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law. Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Attention: To check the authenticity of testing / inspection report & certificate, please contact us at telephone: (86-755) 8307 1443, or email: CN.Doccheck@sgs.com

3rd Building, No. 889, Yishan Road, Xuhui District Shanghai, China 200233
中国·上海·徐汇区宜山路889号3号楼 邮编: 200233

t (86-21) 61402666 f (86-21) 64958763
t (86-21) 61402666 f (86-21) 64958763

www.sgs.com.cn
e sgs.china@sgs.com

Test Report **SL52025233855601TX**
Clause 5.2.5 Microbial Cleanliness
 (EN 14683: 2019 Annex D)

Date: April 20, 2020

Page 3 of 3

	1#	2#	3#	4#	5#
CFU/g	<1	<1	<1	<1	<1

Remark: Performance Requirement: Type I ≤ 30 CFU/g, Type II ≤ 30 CFU/g, Type IIR ≤ 30 CFU/g

The statement of conformity in this test report is only based on measured values by the laboratory and does not take their uncertainties into consideration.

End of Report



Unless otherwise agreed in writing, this document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law. Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Attention: To check the authenticity of testing / inspection report & certificate, please contact us at telephone: (86-755) 8307 1443, or email: CN.Doccheck@sgs.com

SGS-CTC Standards Technical Services (Shanghai) Co., Ltd.
 Testing Center: 3rd Floor, No. 889, Yishan Road, Xuhui District, Shanghai, China

200233
 中国·上海·徐汇区宜山路889号3号楼 邮编: 200233

t (86-21) 61402666 f (86-21) 64958763
 t (86-21) 61402666 f (86-21) 64958763

www.sgs.com.cn
sgs.china@sgs.com

质量管理体系认证证书

证书编号: BMCQ0120088M

兹证明
威海迪尚医疗科技有限公司

统一社会信用代码: 91371000MA3REQTA0Y

注册地址: 山东省威海市火炬高技术产业开发区火炬路 213 号创新创业基地 C 座 406-409 室 邮编: 264200

认证地址: 山东省威海市火炬高技术产业开发区火炬路 213 号创新创业基地 C 座 406-409 室 邮编: 264200

管理体系符合标准:

GB/T 19001-2016 idt ISO 9001:2015

证书产品或服务范围:

医用一次性防护服的生产和销售;隔离衣的生产和销售;一次性使用医用外科口罩的生产和销售;
一次性使用医用口罩的生产和销售;手术衣的生产和销售;医用隔离面罩、一次性使用天然乳胶检
查手套、一次性使用医用 PVC 检查手套、一次性使用医用丁腈检查手套的销售

证书发证日期: 2020 年 06 月 28 日 签发日期: 2020 年 08 月 03 日 证书有效期至: 2023 年 06 月 27 日

注: 证书有效期内, 获证组织须接受发证机构的年度确认

本证书信息可在国家认监委官方网站 (www.cnca.gov.cn)

及本公司网站 (www.benchmarks.com.hk) 上查询

第一次监督:

第二次监督:

第三次监督:




签发人

本臻力行认证(上海)有限公司

地址: 中国(上海)自由贸易试验区富特西一路 355 号 10 层 1009 室 电话: 021-58902328

Quality Management System Certification

Certificate Code: BMCQ0120088M

This is to certify that

Weihai Dishang Medical Technology CO.,LTD.

Organization Code: 91371000MA3REQTA0Y

Registered Address: Room 406-409, Site C, Chuang xin Base, NO.213 Huoju Road, Torch High Technology Development District, Weihai City, Shangdong Province. **Postcode:** 264200

Certification Address: Room 406-409, Site C, Chuang xin Base, NO.213 Huoju Road, Torch High Technology Development District, Weihai City, Shangdong Province. **Postcode:** 264200

This Management System

GB/T 19001-2016 idt ISO 9001:2015

Scope of Certification:

Production and Sales of Disposable Medical Protective Clothing; Production and Marketing of Isolation Gown; Production and Marketing of Disposable Medical Surgical Masks; Production and Marketing of Disposable Medical Masks; Production and Marketing of Surgical Gown/Operating Gown; Sales of Face Shield, Disposable Natural Latex Inspection Gloves, Disposable Medical Vinyl Examination Gloves and Disposable Medical Nitrile Examination Gloves.

Date of issue: 06/28/2020

Date of Surveillance: 08/03/2020

Date of expiry: 06/27/2023

Notes: Within the date of validity, organization must accept the annual confirmation of BENCHMARKS.

1st Periodical Audit:

This certificate information can be querying in CNCA website: www.cnca.gov.cn and BENCHMARKS website:

2nd Periodical Audit:

www.benchmarks.com.hk

3rd Periodical Audit:




Signer

Benchmarks Certification (Shanghai) Co., Ltd.

ADD: Room# 1009, 10F, No. 355, Fute West 1 Road, China (Shanghai) Pilot Free Trade Zone, Shanghai City, China. TEL: 021-58902328

证书附件

质量管理体系认证证书

证书编号：BMCQ0120088M

兹证明

威海迪尚医疗科技有限公司

统一社会信用代码：91371000MA3REQTA0Y

经认证地点如下：

场所名称	场所地址	认证范围
威海迪尚医疗 科技有限公司	威海市火炬高技术产业 开发区火炬路 220-5号 / 264200	隔离衣的生产；一次性使用医用 外科口罩的生产；一次性使用医 用口罩的生产；手术衣的生产

证书发证日期：同证书正本

证书有效期至：同证书正本




签发人

医疗器械质量管理体系认证证书

证书编号: BMCMDQ0120003M

兹证明

威海迪尚医疗科技有限公司

统一社会信用代码: 91371000MA3REQTA0Y

注册地址: 山东省威海市火炬高技术产业开发区火炬路213号创新创业基地C座406-409室 邮编: 264200

认证地址: 山东省威海市火炬高技术产业开发区火炬路213号创新创业基地C座406-409室 邮编: 264200

建立的管理体系符合标准:

YY/T 0287-2017 idt ISO 13485: 2016

证书产品或服务范围:

医用一次性防护服的生产和销售;隔离衣的生产和销售;一次性使用医用外科口罩的生产和销售;一次性使用医用口罩的生产和销售;手术衣的生产和销售;医用隔离面罩、一次性使用天然乳胶检查手套、一次性使用医用PVC检查手套、一次性使用医用丁腈检查手套的销售。

证书发证日期: 2020年06月28日 签发日期: 2020年08月03日 证书有效期至: 2023年06月27日

注: 证书有效期内, 获证组织须接受发证机构的年度确认

本证书信息可在国家认监委官方网站 (www.cnca.gov.cn)

及本公司网站 (www.benchmarks.com.hk) 上查询

第一次监督:

第二次监督:

第三次监督:




签发人

Medical Devices Quality Management Systems Certification

Certificate Code: BMCMDQ0120003M

This is to certify that
Weihai Dishang Medical Technology CO., LTD.
Organization Code: 91371000MA3REQTA0Y

Registered Address: Room 406-409, Site C, Chuang Xin Base, NO.213 Huoju Road, Torch High Technology Development District, Weihai City, Shangdong Province. **Postcode:** 264200

Certification Address: Room 406-409, Site C, Chuang Xin Base, NO.213 Huoju Road, Torch High Technology Development District, Weihai City, Shangdong Province. **Postcode:** 264200
This Management System

certification is valid to the following Standard:

YY/T 0287-2017 idt ISO 13485: 2016

Scope of Certification:

Production and Sales of Disposable Medical Protective Clothing; Production and Marketing of Isolation Gown; Production and Marketing of Disposable Medical Surgical Masks; Production and Marketing of Disposable Medical Masks; Production and Marketing of Surgical Gown/Operating Gown; Sales of Face Shield, Disposable Natural Latex Inspection Gloves, Disposable Medical Vinyl Examination Gloves and Disposable Medical Nitrile Examination Gloves.

Date of issue: 06/28/2020 Date of Surveillance: 08/03/2020 Date of expiry: 06/27/2023

Notes: Within the date of validity, organization must accept the annual confirmation of BENCHMARKS.

1st Periodical Audit:

This certificate information can be querying in CNCA websit: www.cnca.gov.cn and BENCHMARKS website:

2nd Periodical Audit:

www.benchmarks.com.hk.

3rd Periodical Audit:



Jody Gu
Signer

证书附件

医疗器械质量管理体系认证证书

证书编号：BMCMDQ0120003M

兹证明

威海迪尚医疗科技有限公司

统一社会信用代码：91371000MA3REQTA0Y

经认证地点如下：

场所名称	场所地址	认证范围
威海迪尚医疗 科技有限公司	威海市火炬高技术产业 开发区火炬路 220-5号 / 264200	隔离衣的生产；一次性使用医用 外科口罩的生产；一次性使用医 用口罩的生产；手术衣的生产

证书发证日期：同证书正本

证书有效期至：同证书正本



签发人



中国认可
国际互认
检测
TESTING
CNAS L10066

Test Report

Report Number: SSMT-R-2020-01188-02

Sample Name: Disposable medical surgical Mask

Study Title: Skin Irritation Test

Standard: ISO 10993-10:2010

Test facility

Jiangsu Science Standard Medical
Testing Co., Ltd.

C4 Building, No.9 Changyang Road, Wujin
District, Changzhou, Jiangsu, China

Sponsor

Weihai Dishang Medical Technology CO
.,Ltd

Room 406-409, Block C, No.213 Torch Road,
Torch High-tech Industrial Development Zone,
Weihai, Shandong, China

Jiangsu Science Standard Medical Testing Co., Ltd.

C4 Building, No.9 Changyang Road, Wujin District, Changzhou, Jiangsu, China 213161 Tel: (86-519-83587899) Fax: (86-519-83587899) www.jsssmmt.com

Contents

Explanation.....	3
Conclusion.....	4
Study verification and signature.....	5
1.0 Purpose.....	6
2.0 Reference.....	6
3.0 Test and control articles.....	6
4.0 Identification of test system.....	6
5.0 Animal Care and Maintenance.....	7
6.0 Justification of the test system.....	7
7.0 Instruments.....	8
8.0 Experiment design and dose.....	8
9.0 Evaluation criteria.....	10
10.0 Results of the test.....	10
11.0 Deviation statement.....	11
12.0 Record.....	11
13.0 Confidentiality agreement.....	11

Explanation

1. Please apply for rechecking within 15 days of receiving the report if there is any objection.
2. Any erasure or without special testing seal renders the report null and void.
3. The report is only valid when signed by the persons who edited, checked and approved it.
4. The result relate only to the articles tested.
5. The report shall not be reproduced except in full without the written approval of the institute.
6. This experiment was carried out in the sub-site and the address is: No. 68, Yaoluo Road, Wujin District, Changzhou City.

Conclusion

The animal skin irritation test was conducted to assess the potential irritation of the test article or material.

The test sample was extracted with 0.9% sodium chloride injection and sesame oil, respectively. The patches (about 2.5 cm×2.5 cm) which moistened by 0.5 ml extract of test article were directly applied to the rabbit skin for 4 hours. Observation for erythema and edema were conducted at 1 h, 24 h, 48 h and 72 h after removal of the patches.

The primary irritation indexes of the polar and non-polar test group were both calculated to be 0. The test result showed that the extract of the test article did not induce skin irritation in rabbit under the test condition.

Study verification and signature

The study was carried out in accordance with the standard operating procedure. The test process was conducted in compliance with the requirements of CNAS-CL01:2018 (ISO/IEC 17025:2017, IDT) and RB/T 214-2017.

Date Received	2020-05-07
Technical Initiation Date	2020-06-08
Technical Completion Date	2020-06-11
Final Report Completion Date	2020-06-11

Edited by Molly

2020.06.11
Date

Checked by Suri

2020.06.11
Date

Approved by Daisy
Authorized signatory

2020.06.12
Date

Jiangsu Science Standard Medical Testing Co., Ltd.



1.0 Purpose

New Zealand rabbits were used to evaluate the potential of skin irritation of samples under the condition of this test.

2.0 Reference

Biological evaluation of medical devices-Part 10:Tests for irritation and skin sensitization(ISO 10993-10:2010)

Biological evaluation of medical devices-Part 12:Sample preparation and reference materials(ISO 10993-12:2012)

Biological evaluation of medical devices-Part 2: Animal welfare requirements (ISO 10993-2:2006)

3.0 Test and control articles

3.1 Test article (The information about the test article was supplied by the sponsor wherever applicable.)

Test article name: Disposable medical surgical Mask

Sterilization state: Unsterilized

Model: DS01-WHITE/DS02-BLUE

Size: N/S

Lot/ Batch#: N/S

Physical State: Solid

Color: See the photo

Density: N/S

Stability: N/S


Solubility: N/S

Test Article Material: N/S

Packing Material: N/S

Storage Condition: Room Temperature

Sample photograph:

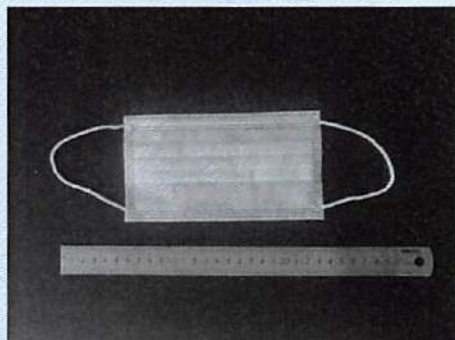

3.2 Control Articles

3.2.1 Polar Negative Control: 0.9% Sodium chloride injection (SC)

Manufacturer: Chenxin Pharmaceutical Co., Ltd.

Size: 250 ml

Physical State: Liquid

Color: Colourless

Lot/ Batch#: 1906112830

Storage Condition: Room Temperature

3.2.2 Non-polar Negative Control: Sesame Oil (SO)

Manufacturer: Ji'an Ivyuanxiangliao. Co., Ltd.

Size: 20 kg

Physical State: Liquid

Color: Pale yellow

Lot/ Batch#: 20190516

Storage Condition: Room Temperature

4.0 Identification of test system

Species: New Zealand white rabbit (single strain)

Number: 6 (3 for polar test group and 3 for non-polar group)

Sex: Female

Weight: Initial body weight not less than 2.0 kg

Health status: Healthy, young adult, nulliparous and not pregnant.

Housing: Animals were housed in groups in cages identified by a card indicating the lab number and test code.

Animal identification: Cage card

The quarantine period: 5 days

5.0 Animal Care and Maintenance

Animal purchase: Provided by Tongxiang Yin Hai Animal Husbandry Professional Cooperative <Permit Code: SCXK (ZHE) 2018-0002>

Bedding: NA

Feed: Rabbit Diet, Beijing Keao Xieli Feed Co., Ltd.

Water: Drinking water met the Standards for Drinking Water Quality (GB 5749-2006)

Cages: Stainless steel cage, Suzhou Fengqiao purification equipment Co., Ltd.

Environment: Temperature 16-26°C, Relative humidity 40%-70%, Lights 12 hours light/dark cycle

Personnel: Associates involved were appropriately qualified and trained

Selection: Only healthy animals were selected

Veterinarian: Vet takes care of the whole course

Ethics: Test methods of operation were reviewed and approved by the Commission on Science Standard animal ethics

There were no known contaminants present in the feed, water expected to interfere with the test data.

6.0 Justification of the test system

6.1 The rabbit is specified as an appropriate animal model for evaluating potential skin irritants by the current

testing standards. Positive control 15% sodium dodecyl sulfate has been substantiated at SSMT with this method. Positive control tests are conducted every six months. The last irritation index of polar test group was 6.0. The last irritation index of non-polar test group was 5.7. The data was from the report SSMT-R-2020-01262-01 (Date: 2020-05-29).

6.2 The test article extract was directly applied to the rabbit skin, which was suggested by the standard.

7.0 Instruments

Digital oscillation incubator (SSMT-300)

Electronic balance (SSMT-075)

Clean bench (SSMT-187)

Straight steel ruler (SSMT-210)

8.0 Experiment design and dose

8.1 Sample preparation

The test article was extracted as Table 1. Extract was checked and used immediately after extraction without the process of filtering, centrifugation, dilution, etc. The pH of the extract was not adjusted prior to testing. The preparation process was aseptic. The control article was prepared under the same condition.


Table 1 Sample Preparation

Aseptic Sampling		Aseptic Agitation Extraction In Inert Container				Final Extract
Sampling Manner	Actually Sampling	Extraction solvent	Extraction ratio	Solvent volume	Condition	Clear or Not
Random sampling	30.0 cm ²	0.9% sodium chloride injection	3 cm ² : 1 ml	10.0 ml	50 °C, 72 h	Clear
	30.0 cm ²	Sesame oil	3 cm ² : 1 ml	10.0 ml	50 °C, 72 h	Clear

8.2 Test method

Use the rabbits with healthy intact skin. Fur is generally clipped on the back of the rabbits 16 h before testing, a sufficient distance on both sides of the spine for application and observation of all test sites (approximately 10 × 15 cm).

Apply 0.5 ml extract of test article or control to 2.5 cm × 2.5 cm absorbent gauze patches, and then apply the patch soaked with the extract of test article or control directly to the skin on each side as shown in Figure 1. And then wrap the application site with a bandage (semi-occlusive) for 4 h. At the end of the contact time, remove the dressings.

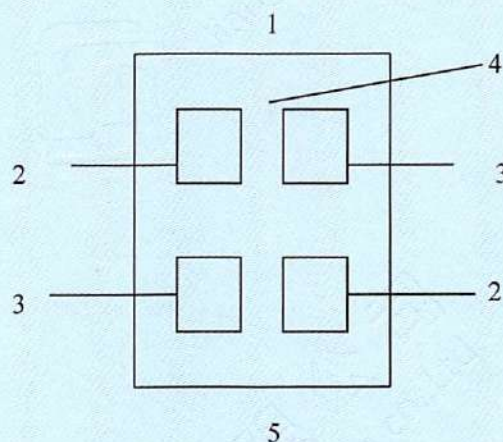

1- Cranial end, 2- Test site, 3- Control site, 4- Clipped dorsal region, 5- Caudal end

Figure1 Location of skin application sites

8.3 Observation of animal

Describe and score the skin reaction for erythema and oedema according to the scoring system given in Table 2 for each application site at each time interval. Record the appearance of each application site at 1 h, 24 h, 48 h and 72 h following removal of the patches.

Table 2 Classification System for Skin Reaction

Reaction	Irritation score
Erythema and Eschar Formation	
No erythema	0
Very slight erythema (barely perceptible)	1
Well-defined erythema	2
Moderate erythema	3
Severe erythema (beet redness) to eschar formation preventing grading of erythema	4
Oedema Formation	
No oedema	0
Very slight oedema (barely perceptible)	1
Well-defined oedema (edges of area well-defined by definite raising)	2
Moderate oedema (raised approximately 1mm)	3
Severe edema (raised more than 1mm and extending beyond exposure area)	4
Maximal possible score for irritation	8

NOTE: Other adverse changes at the skin sites were recorded and are reported.

8.4 Result calculation

Use only 24 h, 48 h and 72 h observations for calculation.

After the 72 h grading, all erythema grades plus oedema grades 24 h, 48 h and 72 h are totalled separately for each test sample and blank for each animal. The primary irritation score for an animal is calculated by dividing the sum of all the scores by 6 (two test/observation sites, three time points).

To obtain the primary irritation index for the test sample add all the primary irritation scores of the individual animals and divide by the number of animals.

Calculate the primary irritation score for the controls and subtract that score from the score using the test material to obtain the primary irritation score.

9.0 Evaluation criteria

The primary irritation index is characterized by number (score) and description (response category) given in Table 3.

Table 3 Primary irritation index categories in a rabbit

Mean score	Response category
0-0.4	Negligible
0.5-1.9	Slight
2.0-4.9	Moderate
5-8	Severe

10.0 Results of the test

According to what observed, the response of skin on testing side did not exceed that on the control side. See Table 4.

Table 4 Dermal observations

Extraction solvent	Rabbit No.	Group		Interval			
				1h	24h	48h	72h
0.9% sodium chloride injection	J1501	Test Article	Erythema	0	0	0	0
			Oedema	0	0	0	0
		Negative Control	Erythema	0	0	0	0
			Oedema	0	0	0	0
	J1502	Test Article	Erythema	0	0	0	0
			Oedema	0	0	0	0
		Negative Control	Erythema	0	0	0	0
			Oedema	0	0	0	0
	J1503	Test Article	Erythema	0	0	0	0
			Oedema	0	0	0	0
		Negative Control	Erythema	0	0	0	0
			Oedema	0	0	0	0
Sesame oil	F1501	Test Article	Erythema	0	0	0	0
			Oedema	0	0	0	0
		Negative Control	Erythema	0	0	0	0
			Oedema	0	0	0	0
	F1502	Test Article	Erythema	0	0	0	0
			Oedema	0	0	0	0
		Negative Control	Erythema	0	0	0	0
			Oedema	0	0	0	0
	F1503	Test Article	Erythema	0	0	0	0
			Oedema	0	0	0	0
		Negative Control	Erythema	0	0	0	0
			Oedema	0	0	0	0

The primary irritation indexes of the polar and non-polar test group were both calculated to be 0. Under the conditions of this study, the extract of the test article did not induce skin irritation.

11.0 Deviation statement

There was no deviation from the standard operating procedure which were judged to have any impact on the validity of the data.

12.0 Record

All the original data and records related to this test and copies of the final report are retained in the archives of Science Standard Medical Testing.

13.0 Confidentiality agreement

Statements of confidentiality were as agreed upon prior to study initiation.





中国认可
国际互认
检测
TESTING
CNAS L10066

Test Report

Report Number: SSMT-R-2020-01188-01

Sample Name: Disposable medical surgical Mask

Study Title: In Vitro Cytotoxicity Test

Standard: ISO 10993-5:2009

Test facility

Jiangsu Science Standard Medical
Testing Co., Ltd.

C4 Building, No.9 Changyang Road, Wujin
District, Changzhou, Jiangsu, China

Sponsor

Weihai Dishang Medical Technology CO.,
Ltd

Room 406-409, Block C, No.213 Torch Road,
Torch High-tech Industrial Development Zone,
Weihai, Shandong China

Jiangsu Science Standard Medical Testing Co., Ltd.

C4 Building, No.9 Changyang Road, Wujin District, Changzhou, Jiangsu, China 213161 Tel: (86-519-83587899) Fax: (86-519-83587899) www.jssmt.com

Contents

Explanation.....	3
Conclusion.....	4
Study verification and signature.....	5
1.0 Purpose.....	6
2.0 Standard.....	6
3.0 Test and control articles.....	6
4.0 Identification of test system.....	7
5.0 Justification of test system.....	7
6.0 Instruments and Reagents.....	7
7.0 Experiment design and dose.....	8
8.0 Evaluation criteria.....	9
9.0 Results of the test.....	9
10.0 Deviation statement.....	10
11.0 Record.....	10
12.0 Confidentiality agreement.....	10

Explanation

1. Please apply for rechecking within 15 days of receiving the report if there is any objection.
2. Any erasure or without special testing seal renders the report null and void.
3. The report is only valid when signed by the persons who edited, checked and approved it.
4. The result relate only to the articles tested.
5. The report shall not be reproduced except in full, without approval of the laboratory.

Conclusion

The study was to investigate the potential cytotoxicity of the test sample. The extract of the test article was added to L-929 cells and then incubated at 37°C in 5% CO₂ for 24 hours. After the incubation, observe the cell morphology. The results were detected with MTT method. The results showed that the cytotoxicity ratio of the 100 % test article extract was 84.5% and the results of control groups showed the test was valid.

Under the conditions of this study, the extract of the test article did not show potential toxicity to L-929 cells.

Study verification and signature

The study was carried out in accordance with the standard operating procedure. The test process was conducted in compliance with the requirements of CNAS-CL01:2018 (IDT ISO/IEC 17025:2017) and RB/T 214-2017.

Date Received	2020-05-07
Technical Initiation Date	2020-05-19
Technical Completion Date	2020-05-21
Final Report Completion Date	2020-05-25

Edited by

Cindy

2020.05.25
Date

Checked by

Bella

2020.05.25
Date

Approved by

Daisy
Authorized Signatory2020.06.12
Date

Jiangsu Science Standard Medical Testing Co., Ltd.



1.0 Purpose

The purpose of the test is to determine the potential cytotoxicity toxicity of a mammalian cell culture (mouse fibroblast L-929 cells) in response to the test article.

2.0 Standard

Biological evaluation of medical devices Part 5: Tests for In Vitro Cytotoxicity (ISO 10993-5:2009)

Biological evaluation of medical devices Part 12: Sample preparation and reference materials (ISO 10993-12:2012)

3.0 Test and control articles

3.1 Test article (The information about the test article was supplied by the sponsor wherever applicable.)

Test article name: Disposable medical surgical Mask

Sterilization state: Not sterilized

Model: DS01-WHITE/DS02-BLUE

Size: N/S

Lot/ Batch#: N/S

Physical State: Solid

Color: See the photo

Density: N/S

Stability: N/S


Solubility: N/S

Test Article Material: N/S

Packing Material: N/S

Storage Condition: Room temperature

Sample photograph:

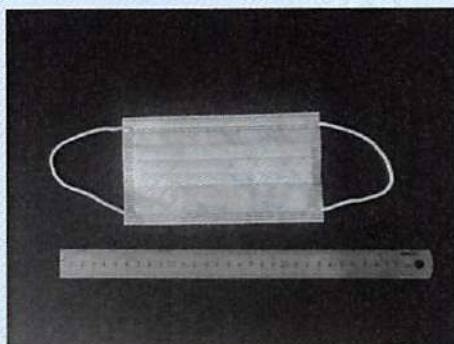

3.2 Control Articles

3.2.1 Negative Control Article Name: High Density Polyethylene

Manufacturer: Jiangsu haiaosihui biotechnology co., LTD.

Size: 1.6 mm thick, 300*300 mm

Lot/ Batch#: M02F017

Physical State: Solid

Color: White

Storage Conditions: Room temperature

3.2.2 Positive Control Article Name: ZDEC

Manufacturer: Tokyo Into Industrial Co., Ltd.

Size: 25 g

Lot/ Batch#: DUDQG-JF

Physical State: Solid

Color: White

Storage Condition: Room temperature

Concentration: 0.1%

3.2.3 Blank Control Name: MEM medium, with addition 10% FBS

Physical State: Liquid

Color: Pink

Storage Condition: 4 °C

4.0 Identification of test system

Mouse fibroblast L-929 cells obtained from ATCC CCL1 (NCTC clone 929).

5.0 Justification of test system

5.1 Historically, mouse fibroblast L-929 cells have been used for cytotoxicity studies because they demonstrate sensitivity to extractable cytotoxic articles.

5.2 The test article was extracted and administered in vitro to mouse fibroblast L-929 cells through a solvent compatible with the test system. This was the optimal route of administration available in this test system as recommended in the standard.

6.0 Instruments and Reagents

6.1 Instruments

CO₂ Incubator (SSMT-279)

Biological microscope (SSMT-278)

Clean bench (SSMT-028)

Bench type low speed centrifuge (SSMT-048)

Vapour-bathing Constant Temperature Vibrator (SSMT-004)

Electronic balance (SSMT-015)

Steel Straight Scale (SSMT-072)

Multiskan Spectrum Microplate Spectrophotometer (SSMT-139)

Mini Vibrator (SSMT-311)

6.2 Reagents

FBS

MEM

Trypsin

Penicillin, Streptomycin sulfate

PBS

MTT (3-(4,5-Dimethylthiazol-2-yl)-2,5-diphenyltetrazolium bromide)

Isopropyl alcohol

7.0 Experiment design and dose

7.1 Sample preparation

Aseptic extracting the test article (test article to volume of vehicle) according to the table below. Sealed and incubated in Vapour-bathing Constant Temperature Vibrator at 37 °C and 60 rpm for 24 hours. After the extraction, check the extraction changes, and immediately use for the experiment, the leach was not filtered, centrifuged or diluted. No pH adjustment.

Table 1 Sample preparation

Sampling		Aseptic Agitation Extraction In Inert Container				Final Extract
Sampling Manner	Actually sampling	Extraction solvent	Extraction ratio	Solvent volume	Condition	Clear or Not
Random sampling	18.0 cm ²	MEM medium (10% FBS)	3 cm ² : 1 ml	6.0 ml	37 °C, 24 h	Clear

The blank control (MEM medium, with addition 10% FBS) and negative/positive controls were prepared in the same condition.

7.2 Test method

Aseptic procedures were used for handling cell cultures.

L-929 cells were cultured in MEM medium (10% FBS, Penicillin 100 U/ml, Streptomycin sulfate 100 µg/ml) at 37°C in a humidified atmosphere of 5% CO₂, then digested by 0.25% trypsin containing EDTA to get single cell suspension. And obtain a 1×10^5 cells/ml suspension by centrifuging (200 g, 3 min) and re-dispersing in MEM medium finally.

The suspended cells were dispensed at 100 µl per well in 96-well plate, and cultured in cell incubator (5% CO₂, 37°C, >90%humidity). Cell morphology was evaluated to verify that the monolayer was satisfactory.

After the cells grew to form a monolayer, original culture medium was discarded. The 96-well plates were then treated with 100 µl of extract of test article (100%, 75%, 50%, 25%), control article, negative article (100%) and positive article (100%) respectively. The 96-well plate was incubated at 37°C in cell incubator of 5% CO₂ for 24 h. Six replicates of each test were tested.

After 24 h incubation, observe the cell morphology first and then discard the culture medium. A 50 µl aliquot of MTT (1 mg/ml) was added to each well and then incubated at 37°C in a humidified atmosphere of 5% CO₂ for 2 hours. The liquid in each well was tipped out and 100 µl isopropanol was added to each well to suspend the cell layer. The microporous plate was vibrated for 10 min and monitored by the optical density at 570 nm on the microplate analyzer.

7.3 Statistical method

Mean±standard deviation ($\bar{x} \pm s$)

Viab. % = $100 \times OD_{570e} / OD_{570b}$

Where: OD_{570e} ——is the mean value of the measured optical density of test sample/negative control/positive control;

OD_{570b} ——is the mean value of the measured optical density of the blanks.

7.4 Observation of the cell morphology

Table 2 Observation of the cell morphology

Grade	Conditions of all cultures
0	Discrete intracytoplasmatic granules, no cell lysis, no reduction of cell growth.
1	Not more than 20 % of the cells are round, loosely attached and without intracytoplasmatic granules, or show changes in morphology; occasional lysed cells are present; only slight growth inhibition observable.
2	Not more than 50 % of the cells are round, devoid of intracytoplasmatic granules, no extensive cell lysis; not more than 50 % growth inhibition observable.
3	Not more than 70 % of the cell layers contain rounded cells or are lysed; cell layers not completely destroyed, but more than 50 % growth inhibition observable.
4	Nearly complete or complete destruction of the cell layers.

8.0 Evaluation criteria

8.1 The 50% extract of the test article should have at least the same or a higher viability than the 100% extract. Otherwise the test should be repeated.

8.2 The lower the Viab.% value, the higher the cytotoxic potential of the test article is.

8.3 If viability is reduced to < 70% of the blank, it has a cytotoxic potential.

8.4 The Viab.% of the 100% extract of the test article is the final result.

9.0 Results of the test

Table 3 Results of the cell vitality

Group	$\bar{x} \pm s$	Viability%	The morphology of the extracted cells was observed under the microscope
Blank control	0.523±0.042	100.0	0
Negative control	0.523±0.010	99.9	0
Positive control	0.013±0.004	2.5	4
100% test article extract	0.442±0.021	84.5	0
75% test article extract	0.469±0.019	89.6	0

50% test article extract	0.494±0.029	94.5	0
25% test article extract	0.519±0.015	99.2	0
Quality check	<p>The mean OD₅₇₀ of blanks is ≥ 0.2.</p> <p>The left (row2) and the right (row11) mean of the blanks do not differ by more than 15 %.</p> <p>The test meets the acceptance criteria.</p>		
Conclusion	Under the conditions of this study, the test article did not show potential toxicity to L-929 cells.		

10.0 Deviation statement

There was no deviation from the approved standard operating procedure which were judged to have any impact on the validity of the data.

11.0 Record

All the original data and records related to this test and copies of the final report are retained in the archives of Science Standard Medical Testing.

12.0 Confidentiality agreement

Statements of confidentiality were as agreed upon prior to study initiation.





中国认可
国际互认
检测
TESTING
CNAS L10066

Test Report

Report Number: SSMT-R-2020-01188-03

Sample Name: Disposable medical surgical Mask

Study Title: Skin Sensitization Test - 0.9% Sodium
Chloride Injection Extract

Standard: ISO 10993-10:2010

Test facility

Jiangsu Science Standard Medical
Testing Co., Ltd.

C4 Building, No.9 Changyang Road, Wujin
District, Changzhou, Jiangsu, China

Sponsor

Weihai Dishang Medical Technology CO
.,Ltd

Room 406-409, Block C, No.213 Torch Road,
Torch High-tech Industrial Development Zone,
Weihai, Shandong, China

Jiangsu Science Standard Medical Testing Co., Ltd.

C4 Building, No.9 Changyang Road, Wujin District, Changzhou, Jiangsu, China 213161 Tel: (86-519-83587899) Fax: (86-519-83587899) www.jssmt.com

Contents

Explanation.....	3
Conclusion.....	4
Study verification and signature.....	5
1.0 Purpose.....	6
2.0 Reference.....	6
3.0 Test and control articles.....	6
4.0 Identification of test system.....	7
5.0 Animal Care and Maintenance.....	7
6.0 Justification of the test system.....	7
7.0 Instruments and reagents.....	8
8.0 Experiment design and dose.....	8
9.0 Evaluation criteria.....	10
10.0 Results of the test.....	10
11.0 Deviation statement.....	10
12.0 Record.....	10
13.0 Confidentiality agreement.....	11

Explanation

1. Please apply for rechecking within 15 days of receiving the report if there is any objection.
2. Any erasure or without special testing seal renders the report null and void.
3. The report is only valid when signed by the persons who edited, checked and approved it.
4. The result relate only to the articles tested.
5. The report shall not be reproduced except in full without the written approval of the institute.
6. This experiment was carried out in the sub-site and the address is: No. 68, Yaoluo Road, Wujin District, Changzhou City.

Conclusion

The extract of the test article was evaluated for its potential to induce skin sensitization in the Guinea Pig Maximization Test.

The test article were extracted with 0.9% sodium chloride injection. The test article extract was intradermally injected and applied topically for induction to ten guinea pigs. Five control animals were treated accordingly but with the solvent alone.

The topical challenge with the extract of test article elicited no skin reaction in the test and the control animals. The skin sensitization rate was determined with 0%.

Study verification and signature

The study was carried out in accordance with the standard operating procedure. The test process was conducted in compliance with the requirements of CNAS-CL01:2018 (ISO/IEC17025:2017, IDT) and RB/T214-2017.

Date Received	2020-05-07
Technical Initiation Date	2020-05-18
Technical Completion Date	2020-06-11
Final Report Completion Date	2020-06-11

Edited by Molly

2020.06.11
Date

Checked by Suri

2020.06.11
Date

Approved by Daisy
Authorized signatory

2020.06.12
Date

Jiangsu Science Standard Medical Testing Co., Ltd.



1.0 Purpose

The test was designed to evaluate the potential of a test article to cause skin sensitization using Guinea Pig Maximization Test.

2.0 Reference

Biological evaluation of medical devices-Part 10:Tests for irritation and skin sensitization(ISO 10993-10:2010)

Biological evaluation of medical devices-Part 12:Sample preparation and reference materials(ISO 10993-12:2012)

Biological evaluation of medical devices-Part 2: Animal welfare requirements (ISO 10993-2:2006)

3.0 Test and control articles

3.1 Test article (The information about the test article was supplied by the sponsor wherever applicable.)

Test article name: Disposable medical surgical Mask

Sterilization state: Unsterilized

Model: DS01-WHITE/DS02-BLUE

Size: N/S

Lot/ Batch#: N/S

Physical State: Solid

Color: See the photo

Density: N/S

Stability: N/S


Solubility: N/S

Test Article Material: N/S

Packing Material: N/S

Storage Condition: Room Temperature

Sample photograph:

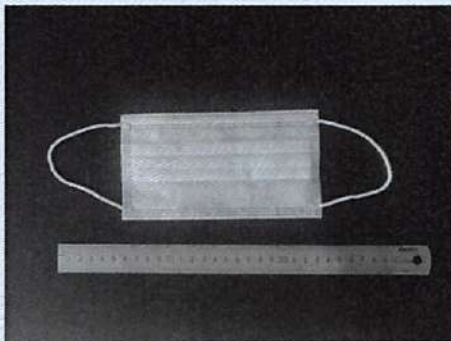

3.2 Control Article

Name: 0.9% Sodium chloride injection (SC)

Manufacturer: Chenxin Pharmaceutical Co., Ltd.

Size: 250 ml

Physical State: Liquid

Color: Colourless

Lot/ Batch#: 1906112830

Storage Condition: Room Temperature

4.0 Identification of test system

Species: Hartley Guinea Pig (*Cavia Porcellus*)

Number: 15 (10 for Test and 5 for Negative Control)

Sex: Males

Health status: Healthy, not previously used in other experimental procedures

Housing: Animals were housed in groups in cages identified by a card indicating the lab number and test code.

Animal identification: Stain with picric acid

The quarantine period: 5 days

5.0 Animal Care and Maintenance

Animal purchase: Provided by Suzhou Genesc Biotechnology Co., Ltd <Permit Code: SCXK (SU) 2015-0002>

Bedding: NA

Feed: Guinea Pig Diet, Beijing Keao Xieli Feed Co., Ltd.

Water: Drinking water met the Standards for Drinking Water Quality GB 5749-2006

Cages: Plastic cage, Suzhou Fengqiao purification equipment Co., Ltd.

Environment: Temperature 18-29°C, Relative humidity 40%-70%, Lights 12 hours light/dark cycle

Personnel: Associates involved were appropriately qualified and trained

Selection: Only healthy, previously unused animals were selected

Veterinarian: Vet takes care of the whole course

Ethics: Test methods of operation were reviewed and approved by the Commission on Science Standard animal ethics

There were no known contaminants present in the feed, water, or bedding expected to interfere with the test data.

6.0 Justification of the test system

6.1 The guinea pig is believed to be the most sensitive animal model for this type of study. The susceptibility of the guinea pig to a known sensitizing agent, dinitrochlorobenzene (DNCB) has been substantiated at SSMT. The skin sensitized positive control test is conducted every six months. The last allergenic rate is 100%. The data was from the report SSMT-R-2020-00198-01 (Date: 2020-04-05).

6.2 The test article was extracted and administered in vivo through a medium compatible with the test system, which is considered as the best route of administration.

7.0 Instruments and reagents

7.1 Instruments

Digital oscillation incubator (SSMT-300)

Steel straight ruler (SSMT-210)

Electronic balance (SSMT-075)

Electronic balance (SSMT-147)

Clean bench (SSMT-187)

7.2 Reagents

Sodium dodecyl sulfate (SDS)

Freund's Adjuvant, Complete liquid

8.0 Experiment design and dose

8.1 Sample preparation

The test article was extracted as Table 1. Extract was checked and used immediately after extraction without the process of filtering, centrifugation, dilution, etc. The pH of the extract was not adjusted prior to testing. The preparation process was aseptic. The control article was prepared under the same condition.

Table 1 Sample Preparation

Aseptic Sampling			Aseptic Agitation Extraction In Inert Container				Final Extract
Sampling Manner	Test phase	Actually Sampling	Extraction solvent	Extraction ratio	Solvent volume	Condition	Clear or Not
Random sampling	Intradermal induction phase I	30.0 cm ²	0.9% sodium chloride injection	3 cm ² : 1 ml	10.0 ml	50 °C, 72 h	Clear
	Topical induction phase II	30.0 cm ²			10.0 ml	50 °C, 72 h	Clear
	Challenge phase	30.0 cm ²			10.0 ml	50 °C, 72 h	Clear

8.2 Test method

8.2.1 Intradermal induction phase I

A pair of 0.1 ml intradermal injections was made for each animal, at the sites (A, B and C) in the clipped intrascapular region as shown in the following Figure 1.

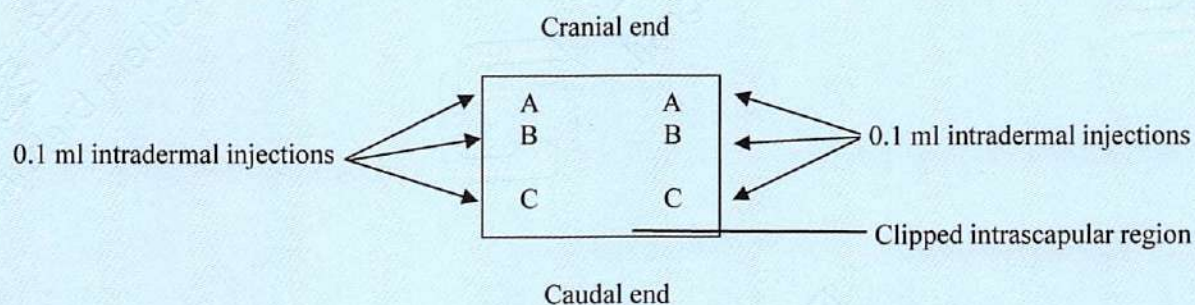


Figure 1 Location of intradermal injection sites

Site A: A 50:50 volume ratio stable emulsion of Freund's complete adjuvant mixed with the solvent.

Site B: The test sample (undiluted extract); inject the control animals with the control articles alone.

Site C: The test sample at the concentration used at site B, emulsified in a 50:50 volume ratio stable emulsion of Freund's complete adjuvant and the solvent (50%); inject the control animals with an emulsion of the blank liquid with adjuvant.

8.2.2 Topical induction phase II

At 6 d after completion of the intradermal induction phase, administer the test sample by topical application to the intrascapular region of each animal, using a patch of area approximately 8 cm² (absorbent gauze) soaked with 0.5 ml extract, so as to cover the intradermal injection sites. Use the concentration selected in the intradermal induction phase for site B. If the maximum concentration that can be achieved in Intradermal induction phase I did not produce irritation, animals were pretreated with 10% sodium dodecyl sulfate 24 hours before the topical induction application. Secure the patches with an occlusive dressing. Remove the dressings and patches after 48 h.

Treat the control animals similarly, using the blank liquid alone.

8.2.3 Challenge phase

At 13 d after completion of the topical induction phase, challenge all test and control animals with the test sample. Administer the test sample and a blank by topical application to left and right abdomen of animals respectively, using absorbent gauze (about 8 cm²) soaked with 0.5ml extracts or solvent control. Secure with an occlusive dressing. Remove the dressings and patches after 24 h.

8.3 Observation of animal

Observe the appearance of the challenge skin sites of the test and control animals 24 h and 48 h after removal of the dressings. Full-spectrum lighting was used to visualize the skin reactions. Describe and grade the skin reactions for erythema and oedema according to the Magnusson and Kligman grading given in Table 2 for each challenge site and at each time interval.

Table 2 Magnusson and Kligman scale

Patch test reaction	Grading scale
No visible change	0
Discrete or patchy erythema	1
Moderate and confluent erythema	2
Intense erythema and/or swelling	3

9.0 Evaluation criteria

Magnusson and Kligman grades of 1 or greater in the test group generally indicate sensitization, provided grades of less than 1 are seen in control animals.

If grades of 1 or greater are noted in control animals, then the reactions of test animals which exceed the most severe reaction in control animals are presumed to be due to sensitization.

If the response is equivocal, rechallenge is recommended to confirm the results from the first challenge.

10.0 Results of the test

The skin response of guinea pigs is shown in Table 3.

Table 3 Guinea pig Sensitization Dermal Reactions

Group	Animal Number	Excitation patch removed 24 h	Excitation patch removed 48 h	Positive rate after challenge phase	Weight range before injection (g)	Weight range after experiment (g)	Abnormal appearance except dermal reactions
Control	J1001	0	0	0%	316.4-379.2	476.5-531.6	None
	J1002	0	0				None
	J1003	0	0				None
	J1004	0	0				None
	J1005	0	0				None
Test	J2001	0	0	0%	309.2-365.2	459.3-539.2	None
	J2002	0	0				None
	J2003	0	0				None
	J2004	0	0				None
	J2005	0	0				None
	J2006	0	0				None
	J2007	0	0				None
	J2008	0	0				None
	J2009	0	0				None
	J2010	0	0				None

Under the conditions of this study, the test article did not show significant evidence of causing skin sensitization in the guinea pig. The skin sensitization rate was determined with 0%.

11.0 Deviation statement

There was no deviation from the standard operating procedure which were judged to have any impact on the validity of the data.

12.0 Record

All the original data and records related to this test and copies of the final report are retained in the archives of Science Standard Medical Testing.

13.0 Confidentiality agreement

Statements of confidentiality were as agreed upon prior to study initiation.

